



**VELOX**<sup>TM</sup>  
**Solvent Safety Cap**




# VELOX<sup>TM</sup>

## Solvent Safety Cap

The extraction of fluids from improperly sealed mobile phase bottles results in the suction of dust, air dirt particles that could be drawn into the capillaries which are not properly in position making their way into highly sensitive analysis systems

Velox Solvent safety caps provides an airtight seal around the bottles with filter function working in dual ways by preventing the dust, contaminant particles from air entering the mobile phase, at the same time blocking the volatile vapour's escaping and hence ensuring constant precise mixing ratios of solvents throughout the time

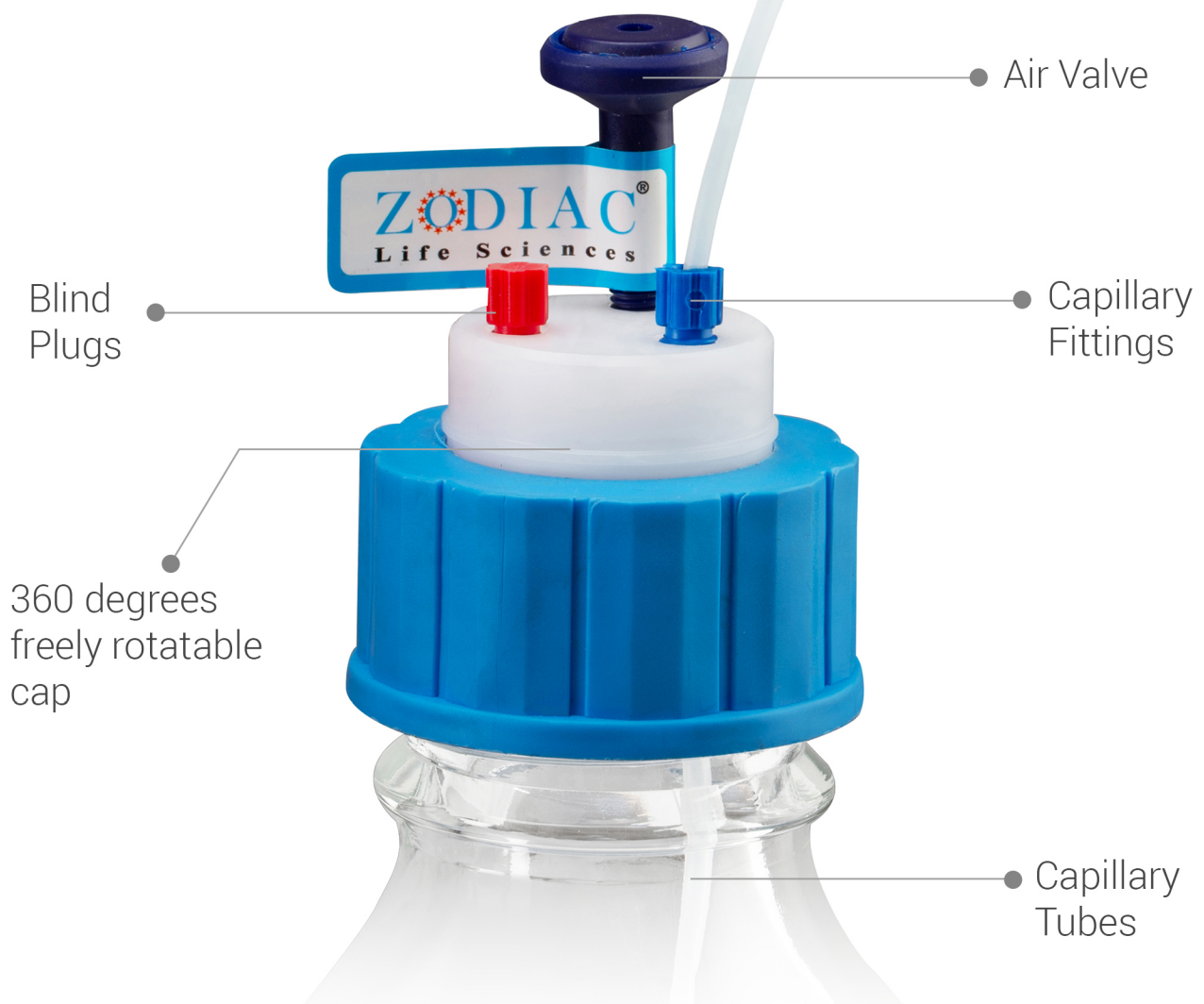

### Advantages

- All harmful vapours are blocked
- Easy container changes with freely rotatable cap
- No accidental intake of air into HPLC system
- With Velox you pass quality and safe inspections with confidence





The new, improved air valve combines valve and filter functions ensuring safe extractions of the solvents so the results of your analysis stay correct

The valve fits your existing Velox Safety Caps without any technical modification

Replace every 6 months for optimal protection as per our filter studies.



- The air valve with the part number is delivered with a change-advice-label.
- Just sign-in the current date at the installation and you are safe in observing the service life.



### ORDERING INFORMATION

Part No.	Description	Thread size	Connector 3.2mm O OD
ZLSIVSC 01	Safety Cap I	GL-45	1X
ZLSIVSC 02	Safety Cap II	GL-45	2X
ZLSIVSC 03	Safety Cap III	GL-45	3X
ZLSIVSC 04	Safety Cap IV	GL-45	4X
ZLSIVSCS 01	Safety Cap Set	GL-45	4X
ZLSAV	Air Valve	GL-45	

